SMOOTHSKIN
pure
adapt

IPL USER GUIDE

ENGLISH

Page number 4

ESPAÑOL

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Welcome to

SMOOTHSKIN PURE ADAPT!

Congratulations! You've opened up a world of unparalleled beauty, smoothness and convenience. Come and join in the conversation at www.smoothskin.com where you can find expert advice and support throughout your treatment.

CONTENTS


1. THE SMOOTHSKIN PURE ADAPT SYSTEM	5	9. THE CONTOUR ADAPTING SKIN SENSORS	26
2. HOW DOES SMOOTHSKIN PURE ADAPT WORK?	6	10. SMOOTHSKIN PURE ADAPT MODES	27
3. IMPORTANT SAFETY INFORMATION	8	11. CLEANING, MAINTENANCE, STORAGE	29
4. POTENTIAL RISKS, SIDE EFFECTS AND SKIN REACTIONS	14	12. TRAVELLING WITH SMOOTHSKIN PURE ADAPT	30
5. WHAT IS THE INTENDED USE OF SMOOTHSKIN PURE ADAPT	16	13. TROUBLESHOOTING	31
6. YOUR SUITABILITY FOR IPL TREATMENT	17	14. WARRANTY	32
7. HOW TO USE SMOOTHSKIN PURE ADAPT	19	15. LABELS AND SYMBOLS	33
8. SKIN TONE SENSOR	23	16. TECHNICAL SPECIFICATIONS	34

1. THE SMOOTHSKIN PURE ADAPT SYSTEM

PACKAGE CONTENTS

Inside the packaging of SmoothSkin Pure Adapt you will find the following:

- Handset and Base Unit
- Power Cable
- User Manual
- Device Storage Bag

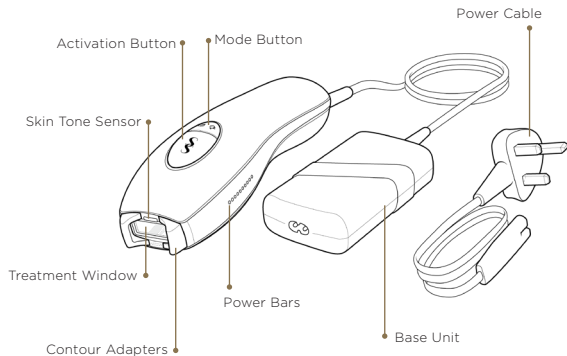
2. HOW DOES SMOOTHSKIN PURE ADAPT WORK?

SmoothSkin Pure Adapt is designed to help break the cycle of hair growth. Light energy is transferred through the skin's surface and is absorbed by melanin present in the hair shaft. The absorbed light energy is converted to heat energy (below the surface of the skin), which reduces the hairs growing back.

Treated hairs naturally fall out over the course of a few days to 1-2 weeks.

Hair growth is cyclical and follows 3 distinct and consecutive phases. IPL treatment is only effective when hairs are in a particular phase (the Anagen - growth phase). Not all hairs are in the same phase at any one time. This is why you should follow the recommended treatment regimen of 12 weekly treatments.

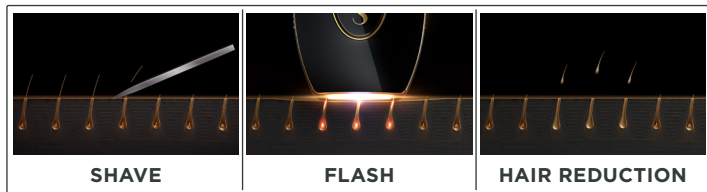WHAT TO EXPECT

Immediately post treatment, you should not see any significant side effects from treatment (Refer to section 4 for more information on side effects).

In the first few weeks following the initial treatments, you will still see some hairs growing. These are likely to be hairs that were missed during treatment i.e. they were lying dormant and not in their growth phase (Anagen) during treatment; when IPL is most effective.

Around 6 weeks into the 12 week treatment programme you should see a reduction in hair growth. However, many hairs may still not have been treated in their growth phase. It is important to continue the weekly treatments.

After the 12 week programme, you should see a significant reduction in hairs within the treated area. Any remaining hairs should be finer and lighter in colour. Continued monthly treatments, or treatments as required, should maintain the reduction in unwanted hair.

3. IMPORTANT SAFETY INFORMATION

CONTRAINDICATIONS

DO NOT use SmoothSkin Pure Adapt if you have very dark skin. The skin tone sensor located near the output window of the device will determine the tone of your skin. If the skin tone sensor determines that your skin is too dark for safe use, operation of SmoothSkin Pure Adapt will be prevented.

DO NOT use SmoothSkin Pure Adapt if your skin has been artificially or naturally tanned or is burned from sun overexposure. Your skin may be extra sensitive following sun exposure and particularly susceptible to burning, blistering, discolouring or scarring if you use SmoothSkin Pure Adapt. Wait 1 week until sunburn or tan disappears before using SmoothSkin Pure Adapt.

DO NOT expose treated areas to the sun. Wait at least 1 week after treatment before exposing treated skin to direct sunlight. Your skin may be extra sensitive following SmoothSkin Pure Adapt treatment and particularly susceptible to sunburn. Apply sunscreen (SPF 15 or greater) to the treated area or cover the treated area with suitable clothing.

DO NOT use SmoothSkin Pure Adapt directly on the nipples, genitals or around the anus. These areas may have a darker skin colour and/or greater hair density. Using the device in these areas may cause discomfort/pain or injure (burn, discolour or scar) your skin.

DO NOT use SmoothSkin Pure Adapt if you have a history of skin cancer or pre-cancerous lesions (e.g. nevi or a large number of moles).

DO NOT use SmoothSkin Pure Adapt if you are pregnant, lactating, as this device has not been tested with these individuals.

DO NOT use SmoothSkin Pure Adapt on dark brown or black spots such as birthmarks, moles or warts in the area you wish to treat. Using SmoothSkin Pure Adapt may injure your skin or make existing conditions worse. You may experience side effects such as burns, blisters and skin colour changes or scarring.

DO NOT use SmoothSkin Pure Adapt if you have any skin condition in the treatment area, including psoriasis, vitiligo, eczema, acne, herpes simplex or active infections or wounds.

DO NOT use on blisters and scars.

DO NOT use SmoothSkin Pure Adapt if you are under the age of 18 as this device has not been tested on these individuals.

WARNINGS

EYE SAFETY

SmoothSkin Pure Adapt emits flashes of intense pulsed light. Direct exposure is potentially harmful to your eyes. Take care to follow the safety precautions below.

WARNING: Possible eye injury (Potentially leading to loss of vision) or skin injury if instructions are not followed. Read and follow the instructions.

The handset can only be activated (flashed) if the skin tone sensor (next to the treatment window) detects a valid skin tone reading and is in full contact with the treatment area. However:

DO NOT look directly into the treatment window of the handset when the system has electrical power and is switched on.

DO NOT attempt to activate (flash) the system towards the eye.

DO NOT treat areas around the eye (eyebrows or eyelashes). Misuse can lead to eye damage, protect eye from exposure. **DO NOT** use the appliance over the eye lids or close to the eye.

ONLY attempt to activate (flash) the system when the front of the handset is in good contact with the area you wish to treat.

WE RECOMMEND you look away from the handset when firing onto your skin.


ELECTRICAL AND SAFETY

As with any electrical device, certain precautions must be taken in order to ensure your safety.

DO NOT try to open SmoothSkin Pure Adapt, as this may expose you to dangerous electrical components.

DO NOT use if SmoothSkin Pure Adapt is damaged; for example, cracked handset, cable damage (internal wires visible), cracked or broken glass on the handset, etc. This may result in injury.

DO NOT use near bathtubs, showers, basins or other vessels containing water. If SmoothSkin Pure Adapt becomes wet, it should not be used.

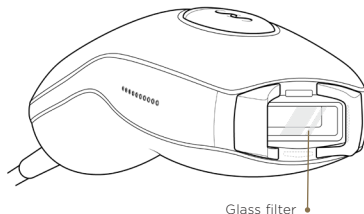WARNING Water and electricity are a dangerous combination!

To avoid risk of electrocution:

- do not use this equipment in wet surroundings (e.g. in the bathroom or near a shower or swimming pool)
- do not let water run into the appliance.

DO NOT use if SmoothSkin Pure Adapt becomes uncomfortably hot to touch. This may indicate that the unit is damaged.

DO NOT use a power supply or other component, other than those provided for use with SmoothSkin Pure Adapt. This may damage the unit and cause it to stop working.

DO NOT attempt to activate (flash) SmoothSkin Pure Adapt onto any surface other than skin.

WARNING no modification of this equipment is allowed. Pure Adapt cannot be serviced

DO NOT position SmoothSkin Pure Adapt so that it is difficult to disconnect from the mains supply.

WARNING: If the wall socket used to power the equipment has poor connections, the plug of the equipment may become hot. Make sure you plug the equipment into a properly installed wall socket to avoid the risk of fire and burns.

GENERAL SAFETY

DO NOT use if you have had a skin peel treatment on the area you wish to treat. Please wait 30 days after a skin peel before treating with SmoothSkin Pure Adapt.

Treating tanned skin with Intense Pulsed Light (IPL) may result in permanent Hyperpigmentation

(Darkening) or Hypopigmentation (lightening) of the treated skin.

DO NOT use on areas where you have tattoos or permanent make up.

DO NOT use If you have a known sensitivity to sunlight (photosensitivity) or are taking medication that makes the skin more sensitive and causes photosensitivity, e.g. Retin A, Accutane and/or other topical retinoids. Always check the instruction leaflet that comes with your medicinal product, to see if photosensitivity is listed as a side effect. Use of SmoothSkin Pure Adapt on photosensitive skin could result in skin damage such as swelling or blisters. Seek medical advice prior to application of the appliance as use may harm the skin.

This appliance is not intended for the treatment of medical conditions. Medical advice should be sought to address, for example, moles, skin rash, itchy skin, skin fungus or infection, skin bumps, or skin tags.

SmoothSkin Pure Adapt should not be used over any fillers or botox.

WARNING: SmoothSkin Pure Adapt should be kept out of the reach of children under 18 years of age. Children should not use, play with, maintain or clean the device. There are a number of potential risks - exposure to the light output causing burns or eye damage, electrocution, strangulation from the cables.

SmoothSkin Pure Adapt should be kept out of reach of pets and pests which could damage the device. Potential risks include reduced function and electrocution.

WARNING: Do not leave the equipment unattended when it is switched on to avoid the risk of fire or burns.

SmoothSkin Pure Adapt is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge. Unless they have been given supervision or instruction concerning use of the equipment to avoid the risk of fire or burns.

SmoothSkin Pure Adapt glass filter is a critical component that will get hot during use.

DO NOT touch the filter during or directly after use. Regularly inspect the filter for damage.

DO NOT use SmoothSkin Pure Adapt if the filter is cracked or missing.

Ensure that SmoothSkin Pure Adapt treatment window is clean and free from debris. Any debris on or around the treatment window during operation may result in temporary skin changes such as redness or swelling.

WARNING: Damage of the skin may occur after prolonged or repeated surface application on one site. Avoid over-usage which included excessive passes, stacking of pulses and/or increased frequency of use.

WARNING: Do not override the safety mechanisms inherent to the device.

DO NOT exceed the usage specified in the treatment regime.

DO NOT use on areas that have recently received local anaesthesia.

PRECAUTIONS

Hair removal by lasers or intense pulse light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern, and South Asian heritage treated on the face and neck.

DO NOT use SmoothSkin Pure Adapt on any area where you may want hair to grow back.

DO NOT use SmoothSkin Pure Adapt device on man's face, jaw, or neck. The hair in a man's beard area is too dense and using SmoothSkin Pure Adapt device in these areas may result in skin injury.

SmoothSkin Pure Adapt is not recommended for use on red, grey, white or blonde hair. SmoothSkin Pure Adapt is not effective on these hair colours.

CAUTION: Do not subject the equipment to heavy shocks to avoid risk of damage to the lamp

CAUTION: Always unplug the appliance after use and in case of power failure to avoid risk of damage to the equipment.

4. POTENTIAL RISKS, SIDE EFFECTS AND SKIN REACTIONS

Use of SmoothSkin Pure Adapt could result in side effects.

The table below lists the known skin reactions that may occur following treatment with the SmoothSkin Pure Adapt.

POSSIBLE SIDE EFFECT	HOW TO ASSESS & REACT
Mild pain / discomfort in the area being treated	This is expected and is normal for all IPL treatments. You can keep on using the device as instructed, and the pain should diminish with continued use.
Warm feeling or tingling sensations during treatment, which typically disappear after a few seconds to a minute and decrease with continued use.	This is expected and is normal for all IPL treatments. You can keep on using the device as instructed.
Itchiness in the treatment area.	This is quite common for IPL treatments, and should subside after a short period. You can keep on using the device as instructed. Do not scratch the area.
Skin redness during or after treatment which disappears within minutes.	This is expected and is normal for all IPL treatments. You can keep on using the device as instructed once skin redness has disappeared.
Skin redness which does not disappear within minutes and lasts 24-48 hours.	Stop using the device immediately and consult your physician before using it again.
Pain or discomfort that is intense during treatment or persists after a treatment.	Stop using the device immediately and consult your physician before using it again.

IN VERY RARE CASES	HOW TO ASSESS AND REACT
Swelling and redness that does not disappear within two to three days.	Stop using the device immediately and consult your physician before using it again.
Temporary changes in skin colour (lightening or darkening).	If your skin color changes, stop using the device immediately and consult your physician.
Blistering, scarring or burning of the skin.	Stop using the device immediately and consult your physician before using it again. Apply a cold pack to the area. Treat with antiseptic or burn cream.

5. WHAT IS THE INTENDED USE OF SMOOTHSKIN PURE ADAPT

INTENDED USE

SmoothSkin Pure Adapt is intended for the permanent reduction in unwanted hair.

INDICATIONS FOR USE

The SmoothSkin Pure Adapt Hair Removal Device is indicated for the removal of unwanted hair. The SmoothSkin Pure Adapt is also indicated for the permanent reduction in hair regrowth, defined as the long-term, stable reduction in the number of hairs regrowing when measured at 6, 9 and 12 months after the completion of a treatment regime.

6. YOUR SUITABILITY FOR IPL TREATMENT

BODY AREAS

SmoothSkin Pure Adapt is suitable for use on body hair and female facial hair below the cheek line.

SKIN COLOUR

SmoothSkin Pure Adapt is suitable for use on light, medium and dark skin tones, up to and including skin tone 5.

SmoothSkin Pure Adapt incorporates a skin tone sensor that measures the colour of

your skin and will prevent use if your skin tone is too dark.

Dark Skin tones are rich in melanin which could absorb too much energy and potentially lead to skin tone damage.

SKINTONE

HAIR COLOUR

SmoothSkin Pure Adapt is suitable for use on naturally black or brown hair. It may not be as effective on white, grey, blonde or red hair.

HAIR COLOUR

7. HOW TO USE SMOOTHSKIN PURE ADAPT

WARNING: Follow each step in the treatment process ensuring you complete all aspects of each step before progressing to the next.

STEP 1: PREPARE THE AREA TO BE TREATED

Remove all visible hair by shaving. Make sure there is no hair remaining above the surface of the skin as this may cause injury. This also prevents debris from covering the front of your device.


Cleanse the treatment area and pat dry. **DON'T** use any gels or creams with this device!

STEP 2: DO A PATCH TEST

Before your first treatment on each new body area, we RECOMMEND you test your skin in that area for a reaction to SmoothSkin Pure Adapt.

The patch test area should be approximately 3cm x 2cm in size (equivalent to 2 flashes applied to the skin side by side). Follow Step 3 to treat the area.

Wait 24 hours following the patch test, to ensure your skin is suitable for treatment and there is no adverse reaction to the light energy. If there is no reaction after 24 hours, you may treat the area around the patch test. The 'patch tested area' should NOT be re-treated for at least 1 week.

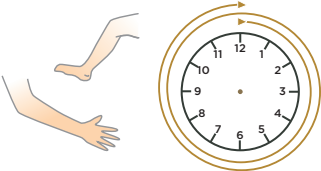An illustration showing two human arms on the left. In the center is a circular clock face with numbers 1 through 12. Two concentric orange circles with arrows indicate a 24-hour cycle. To the right of the clock, the text reads: **24 HOURS**
Observe for adverse events.
Always perform a patch test on each new treatment area

24 HOURS
Observe for adverse events.
Always perform a patch test on each new treatment area

STEP 3: TREATING WITH SMOOTHSKIN PURE ADAPT


- a) Plug the power supply into an electrical outlet. The fans will start up and the skin tone sensors will light up blue, indicating the device is switched on.

POWER MODE Timeout

When the unit is left “inactive” i.e. no buttons are pushed for 10 minutes then the unit will revert to STANDBY mode, and the skin sensor lights will switch off. To return to READY mode for treatment, press the Power Mode button.


- b) Press the handset firmly against the area to be treated ensuring the Skin Tone Sensor and the Skin Contact Sensors are in full contact with skin.
If the skin tone sensor does not detect a valid skin tone and / or the skin contact sensors are not in full contact with skin, then the handset cannot be activated.

When the skin tone sensor measures a valid skin tone and the skin contact sensors are in full contact with skin the Power Bars will illuminate WHITE to indicate the device is ready to flash.

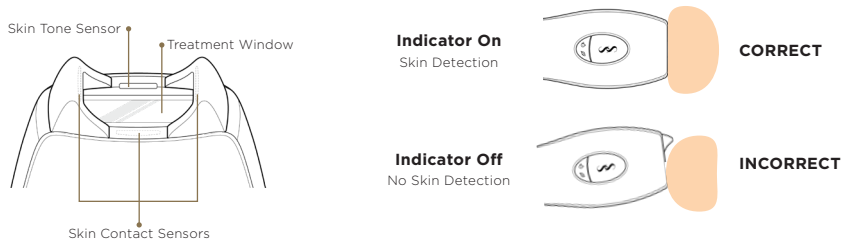- c) Press the activation button. The device will flash, and a warm sensation should be felt.
- d) Move the handset to the next area to be treated and repeat. The device needs to recharge between flashes.

WARNING: Make sure not to cover the vents on the device whilst treating, as this can lead to the device overheating.

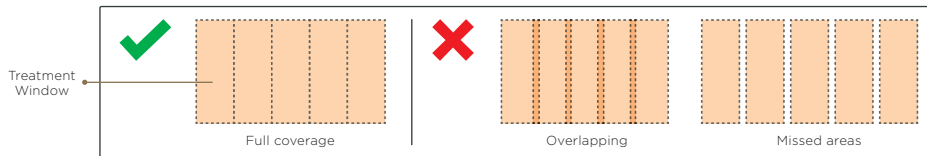TREATMENT METHOD - STAMP OR GLIDE

Stamp Method

This method is recommended for small areas such as the underarm and bikini line. Place the handset on the skin, press and release the activation button. Remove from the skin before carefully re-applying to the next treatment area

Glide Method

This method is recommended for larger areas such as legs. Place the handset on the skin, hold down the activation button (press down continuously), then glide the handset along the skin between flashes. Try to move the handset at a steady speed that gives full coverage without overlap or missed areas.

STEP 4: YOUR SMOOTHSKIN PURE ADAPT TREATMENT SCHEDULE

SmoothSkin Pure Adapt should be used once a week until the desired results are achieved.

STEP 5: SWITCHING OFF SMOOTHSKIN PURE ADAPT

When the treatment is finished, SmoothSkin Pure Adapt should be unplugged from the electrical outlet.

Refer to section 11 for details of how best to look after your SmoothSkin Pure Adapt.

8. SKIN TONE SENSOR

The skin tone sensor measures the colour of your skin and sets the correct level of output power needed. If your skin is too dark for treatment with SmoothSkin Pure Adapt, the skin tone sensor will detect this and prevent the unit from operating. RED lights will be shown in the power bar if your skin is too dark.


The status of the skin tone sensor can be seen on the power bars located on both sides of the handset. If a valid skin tone is detected, the power bars will illuminate WHITE. The number of lights illuminated will be dependent upon your exact skin tone.

The amount of light energy delivered for each of the device settings when in Power Mode (depending on skin colour) is shown in the table below.

This full contact with your skin ensures you get the best results, as all the light energy emitted from the device reaches the treatment area.


NUMBER OF POWER BAR LIGHTS ILLUMINATED	MAXIMUM ENERGY (J/CM²)
10	6.00
9	5.78
8	5.56
7	5.33
6	5.11
5	4.89
4	4.67
3	4.44
2	4.22
1	4.00

POWER BARS

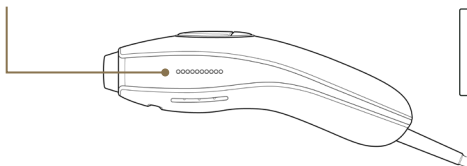
Less lights = Darker skin and lower energy setting
More lights = Lighter skin and higher energy setting


Skin too dark
No treatment possible
1 RED light illuminated

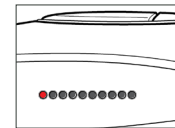Dark skin
Lowest energy setting
1 light illuminated

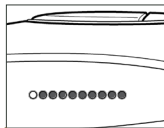Light skin
Highest energy setting
10 lights illuminated

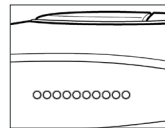Not Safe

Safe and Effective


9. THE CONTOUR ADAPTING SKIN SENSORS

SmoothSkin Pure Adapt has new *Adapt Plus Technology* that adapt to the shape of the treatment area and check for skin contact before each flash, making sure the device is safe to operate! This makes the device easy to use and compliant with the most stringent safety regulations.

The *Adapt Plus Technology* also ensure you get the best results by preventing light leakage and directing the light into your skin.

If the contour adapters are not in full contact with the treatment area, the device will not flash. Try adjusting the positioning of the device to ensure full contact is achieved, particularly on more awkward areas like arms and shins.

Do not try to remove or apply excessive force to the contour adapters as this will damage the device.

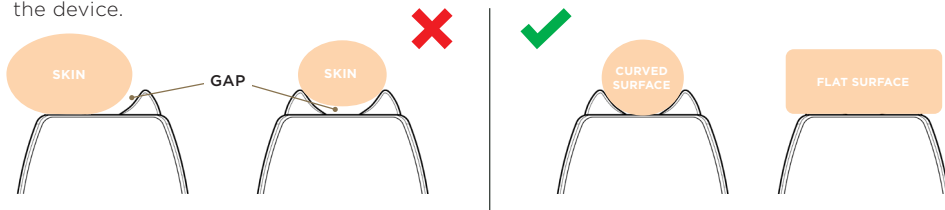10. SMOOTHSKIN PURE ADAPT MODES

You can treat your body using any one of the three SmoothSkin Pure Adapt treatment modes:

Power Mode
No Indicator

Gentle Mode
Indicator On

Speed Mode
Indicator On

Power Mode is the highest power setting of the device. Power Mode will deliver the best results, as it is using the maximum power output possible for your unique skin tone. This mode works particularly well on stubborn areas where the hair tends to be denser, such as underarms and bikini line.

Gentle Mode can be used if you find Power Mode uncomfortable or a little too painful. Operating the device in Gentle Mode will reduce the power output, and you will notice less Power Bar indicators illuminated when applied to your skin.

Speed Mode is the fastest mode, and therefore the lowest power mode. Speed Mode is best for treating larger areas where the hair is less dense (such as arms or legs), and where you want to move along the treatment area quickly. It can also be used as a maintenance mode, to be used to keep hair regrowth at bay after your initial treatment regime.

Power Mode is the default setting when your device is switched on. You can toggle between the three modes pressing the Mode Button underneath the main Activation Button.

When in any of these modes, you can also choose whether to use Glide Method or Stamp Method.

Stamp Method is recommended for small areas such as the underarm and bikini line. Place the handset on the skin, press & release the activation button. Remove from the skin before carefully re-applying to the next treatment area.

Glide Method is recommended for larger areas such as legs. Place the handset on the skin, hold down the activation button (press down continuously), then glide the handset along the skin between flashes. Try to move the handset at a steady speed that gives full coverage without overlap or missed areas.

11. CLEANING, MAINTENANCE, STORAGE

After treatment, always switch off SmoothSkin Pure Adapt by unplugging the electrical outlet.


After use, the handset, skin tone sensor and treatment window should be inspected for damage and wiped down with a dry, lint free cloth.

WARNING: The glass filter can get hot during use. **DO NOT** touch or clean the filter for at least 5 minutes after use to allow it to cool down.

NEVER use water or other cleaning fluids, as these can damage the device leading to a potential safety hazard.

To ensure maximum performance of your SmoothSkin Pure Adapt, it is particularly important to keep the front of the device, Contour Adapters and reflector free from dirt, hair or other debris.

If necessary the reflector and skin tone sensor can be carefully cleaned by using a slightly damp cotton bud.

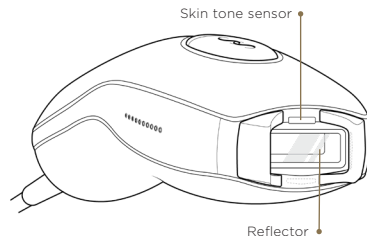Store SmoothSkin Pure Adapt in a cool, dry place. Make sure that the treatment window and skin tone sensor on the handset are protected from damage.

Regularly check the device (including cords) for visible signs of damage. In case of damage or cracks, stop using the device - and visit the website for support details

To prevent strangulation or tripping hazard, always unplug the device from power socket when not in use.

12. TRAVELLING WITH SMOOTHSKIN PURE ADAPT

The power supply unit of SmoothSkin Pure Adapt is universal and will work with all commonly encountered mains supplies. Providing you have a suitable adaptor for the mains power cable plug, SmoothSkin Pure Adapt will operate as normal.

13. TROUBLESHOOTING

PROBLEM	SOLUTION
The power bar indicators on the side of the handset do not light up when placed against your skin.	Make sure the power supply is plugged in. Make sure the handset is in READY mode (press the Standby/Ready button once). Ensure the skin sensors are pressed against the skin. Ensure the skin tone sensor is pressed against the skin.
Whilst pressed against your skin, the first power bar indicator on the side of the handset is illuminated RED, and SmoothSkin Pure Adapt will not work.	The area you are trying to treat is too dark for treatment. SmoothSkin Pure Adapt includes an intelligent skin tone sensor that automatically detects your skin tone and selects the optimal treatment setting. Due to the nature of Intense Pulsed Light technology, SmoothSkin Pure Adapt is not for everyone as it cannot be used on skin tone Type 6 - very dark skin. (See Section 6).
The power bar indicators on the side of the handset are on or are flashing in an unfamiliar pattern including a red LED indicator light.	SmoothSkin Pure Adapt will display any error messages through combinations of lights on the power bars on the sides of the handset. This indicates that there is a fault with the device. Stop using the device and refer to the website.
SmoothSkin Pure Adapt handset, treatment window, power supply or cable are damaged, broken, cracked or appear defective.	DO NOT USE. If you are in any doubt about the safety of SmoothSkin Pure Adapt or suspect it is damaged in any way, it must not be used, and you should refer to the website.

14. WARRANTY

We grant a 2 year warranty on the product starting on the date of purchase. Within the warranty period, we will check faults in materials or workmanship free of charge and replace the complete device, at our discretion.

This warranty extends to every country where this device is supplied by CyDen or its appointed distributor.

This warranty does not cover damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The warranty becomes void if repairs are undertaken by unauthorised persons and if original SmoothSkin Pure Adapt parts are not used. To obtain service within the warranty period, please visit our website **www.smoothskin.com**

15. LABELS AND SYMBOLS

 Type BF Applied Part	 UK Conformity Marking
 Class II Equipment	 European Conformity Marking
 Follow Instructions for use	 FCC EMC compliance for North America
 Keep Dry	 Australia Regulatory Compliance Mark
 WEEE waste disposal mark of the European Union	 Ukraine Compliance Marking
 Warning Optical Radiation from device May Cause Eye Injury	 China RoHS Environmental Protection Use Period
 MET Mark United States & Canada	

FCC Statement:

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

16. TECHNICAL SPECIFICATIONS

SmoothSkin Pure Adapt is a filtered broadband Intense Pulsed Light system with the following technical specifications:

- **Repetition Rate:** Manually pulsed every 0.45 to 1 second (dependent on treatment setting), Continuous operation
- **Pulse:** Partial Free discharge FWHM = 2.5 ± 0.5 ms
- **Max Optical Output:** 3 - 6J/cm²
- **Wavelength:** 510 - 1100nm
- **Power Line Input:** 100-240V~, 50/60Hz, 1.7-0.91A
- **Operating temperature:** Between 5°C and 40°C
- **Operating humidity:** Up to 93%R.H. non condensing
- **Operating Pressures:** 700hPa to 1060hPa
- **Treatment area (spot size):** 3 cm² (30mm x 10mm)

The area directly around the treatment window is classified as the Applied Part, and can reach a maximum temperature of 60°C.

Ensure SmoothSkin Pure Adapt is used as directed in Section 7.

STORAGE & TRANSPORT CONDITIONS:

- **Temperature:** -25°C to +70°C
- **Humidity:** Up to 93%R.H. non condensing
- **Pressures:** 500hPa to 1060hPa

If the device has been stored at temperatures outside of the stated operating range, then please allow at least 30 minutes before use to allow the device to acclimatise.

SmoothSkin Pure Adapt has unlimited flashes and is intended for a single user.

Minimum expected service life is 10 years when following the recommended regimen.

CE marking certifies that this appliance conforms to the following EU Directives:

- Low Voltage Directive 2014/35/EU.
- Electromagnetic Compatibility Directive 2014/30/EU.
- Energy related Product Directive 2009 / 125 / EC.


SmoothSkin Pure Adapt complies with the latest EMC Standards and Amendments applicable to Household and Medical Devices. These are: EN 55014-1:2021, EN 55014-2:2021, EN 61000-3-2:2019+A1:2021, EN 61000-3-3:2013+AMD1:2019+AMD2:2021, IEC 60601-1-2:2014 (Ed 4.0), FCC CFR 47 Part 15.107 & 15.109, ICES-003 Issue 7.

SmoothSkin Pure Adapt complies with the latest Safety Standards and Amendments applicable to Household and Medical Devices.

These are: IEC 60335-1:2010 + Ams, EN 62233, IEC 60335-2-113:2016 IEC 60601-1:2005 (Ed 3.1) +A2:2012, IEC 60601-2-57:2011, IEC 60601-2-83:2019, IEC 60601-1-6:2010 + Amd1:2013, IEC 60601-1-11:2015, IEC 62366-1:2015, IEC 62304:2006 + Ams, CAN/CSA - C22.2 No 60601-1:2008 and 2014, ANSI/AAMI ES 60601-1:2005 + Ams.

SmoothSkin Pure Adapt is classified as Risk Group Exempt for all Hazards.

For further details contact CyDen Ltd.

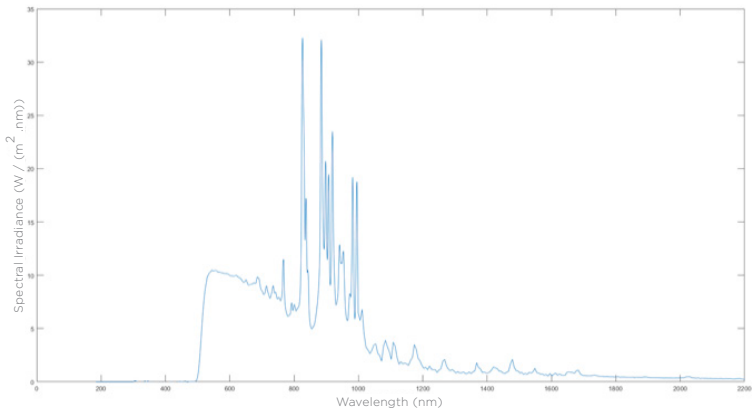DISPOSAL (END OF LIFE)

In order to minimize hazards to health and the environment and ensure that materials can be recycled, this product should be disposed of at a separate collection facility for waste electrical and electronic equipment.

The WEEE symbol marked on the product is there to remind you.

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Póngase en contacto con CyDen Ltd para mayor información.
CyDen Limited, Block A, Bay Studios Business Park, Fabian Way, Swansea SA1 8QB UK
Cyden Ltd. Factory,